



扫码关注公众号



杭州华安生物技术有限公司

☎ 0571-88062880

🖨 0571-88060400

✉ sales@huabio.cn

📍 杭州市钱塘区白杨街道6号大街452号高科技企业孵化园区2号楼16层



ELISA 试剂盒与免疫学检测 ELISA Kits & Immunoassays

- Quality ELISA 试剂盒
 - DTSet ELISA 试剂盒
 - 配对抗体
- 让研究变得简单



HUABIO

为满足不同阶段和精度的蛋白定量需求，华安生物（HUABIO）提供两大ELISA解决方案及其配套免疫学检测试剂：Quality ELISA系列作为开盒即用的标准化工具，主打高效、灵敏与可靠性，适用于大多数定量检测场景；DTset ELISA系列则为需要高度灵活性、大规模筛查或定制化服务的研究者设计。

目录

CONTENTS

- | | |
|----------------------|-------|
| 1. Quality ELISA 试剂盒 | 01-10 |
| 2. DTSet ELISA 试剂盒 | 11-13 |
| 3. 配对抗体 | 14-27 |

完全即用型： Quality ELISA 试剂盒

华安生物（HUABIO）基于双抗夹心法开发了一系列性能稳定，品质卓越的完全即用型Quality ELISA 试剂盒产品。为满足客户对产品质量的深层次需求，在生产、研发和质控环节，华安生物（HUABIO）严格筛选兔单克隆抗体对，并经过多种生物样本验证，确保Quality ELISA试剂盒的最优性能。



1

产品优势

开盒即用

实验流程简化高效，抗体孵育时间低至30min，一天可以做多轮实验

高灵敏度

准确检测pg级别蛋白，低至0.03pg/ml，适用于低丰度蛋白检测

高精密度

CV值低于8%，保证检测结果的一致性

高稳定性

2-8℃保存，有效期高达18个月

CE认证


进口品质，畅销国外，深受国内、外客户认可

实验原理


基于经典的抗原 - 抗体特异性结合与酶联信号放大原理，Quality ELISA系列提供即用型、高性能的蛋白定量解决方案。为适配不同样本类型与实验通量需求，华安生物（HUABIO）开发了“常规系列”与“一步法”两种形式产品。“常规系列”提供经典的灵敏度和宽动态范围；“一步法”则进一步简化操作步骤，在保证高特异性的前提下，为大规模样本筛查提供极致效率。

Quality ELISA 试剂盒的实验流程图


常规系列 抗体孵育时间 90min

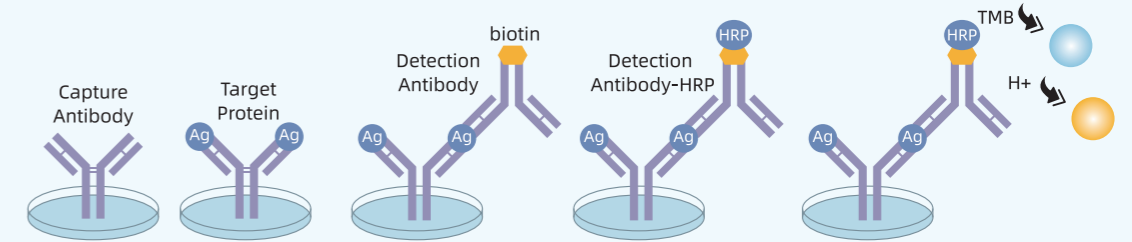

“一步法”系列 抗体孵育时间 30min

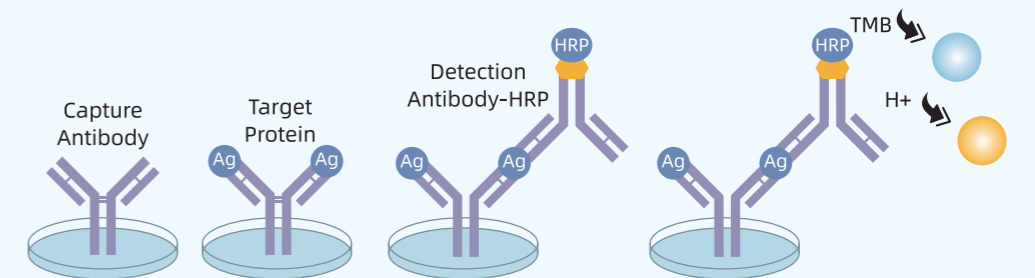

应用领域广，适用于热门领域相关靶点的定性、定量检测

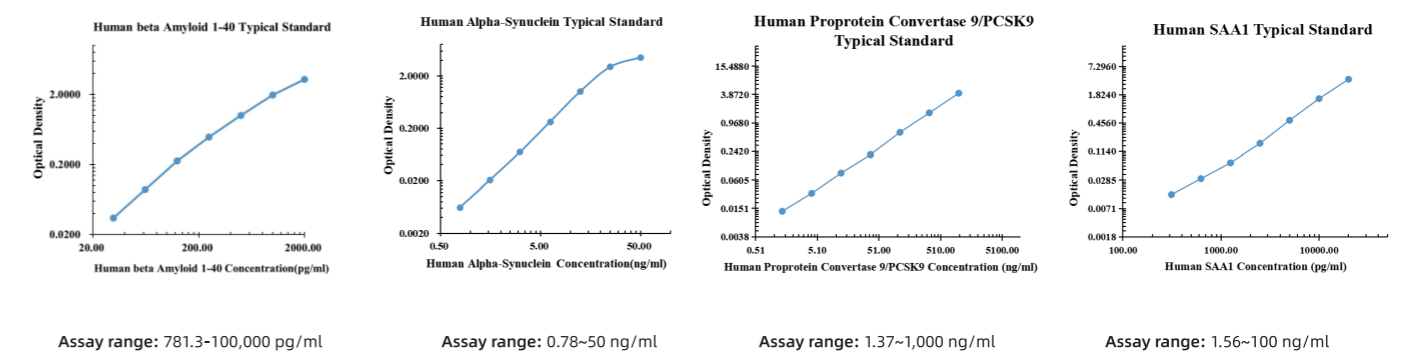



高灵敏度，适用于pg/mL级别蛋白的检测

靶点名称	华安生物 & 攸克生物		竞品 S
	Cat Number	Sensitivity	Sensitivity
Human IL-6	EH0001	0.26 pg/ml	0.7 pg/ml
Human TNF-alpha	EH0002	0.21 pg/ml	6.23 pg/ml
Human IFN-gamma	EH0039	0.36 pg/ml	5.69 pg/ml
Human IL-17A	EH0004	0.26 pg/ml	15 pg/ml
Human IL-4	EH0007	0.03 pg/ml	10 pg/ml

在相同实验操作条件下，对人源IL-6, TNF-alpha等ELISA 试剂盒和竞品S相应靶标进行最低检测限评价实验。结果表明，华安生物（HUABIO）的ELISA 试剂盒灵敏度显著优于竞品，高达竞品的300倍。

精密度优秀，CV值低于8%，保证检测值的一致性

Human IL-6	华安生物 & 攸克生物			竞品 S		
	S1	S2	S3	S1	S2	S3
Sample Number	22	22	22	20	20	20
Average (pg/ml)	24	115.4	263.2	15.8	95.6	179
Standard Deviation	1.1	4.5	9.3	0.7	3	3.1
Coefficient of variation(%)	5.6	3.9	3.6	4.4	3.1	1.7

*精密度评价标准：及格，CV≤15%；良好，CV≤10%；优秀，CV≤8%；卓越，CV≤5%

在相同实验操作条件下，使用人源IL-6 ELISA 试剂盒和竞品S检测3个批次样本中 IL-6 的含量，并进行精密度统计分析。结果表明，华安生物（HUABIO）的ELISA 试剂盒CV值低于8%，严格控制批内/批间差，与竞品S精密度相当。

准确度优秀，线性范围广，确保实验结果的准确可靠

靶点名称	华安生物 & 攸克生物		竞品 S
	Cat Number	Recovery (%)	Recovery (%)
Human IL-6	EH0001	84-115	86-99
Human TNF-alpha	EH0002	95-105	95-108
Human IFN-gamma	EH0039	89-113	90-109
Human IL-17A	EH0004	80-112	90-105
Human IL-4	EH0007	88-105	82-111

*准确度评价标准：及格，回收率在50%~150%；优秀，回收率在80%~120%

在相同实验操作条件下，对人源IL-6, TNF-alpha等ELISA 试剂盒和竞品S相应靶标进行准确度评价分析。结果表明，华安生物（HUABIO）的ELISA 试剂盒回收率严格控制在80%-120%之间，较为优秀，与竞品准确度相当。

客户引用案例（特色文献精选）

Pathogen-derived glyoxylate inhibits Tet2 DNA dioxygenase to facilitate bacterial persister formation

期刊: Cell Metab IF: 27.7 DOI: 2025 May 6;37(5):1137-1151.e5.
产品名称: Mouse CXCL9 ELISA Kit | Catalog# EM0015 客户单位: 复旦大学

Deubiquitinase-dependent transcriptional silencing controls inflammation

期刊: Cell Res IF: 25.9 DOI: 2025 Sep;35(9):675-686.
产品名称: Mouse CXCL2 ELISA Kit | Catalog# EM0018; Mouse IL-1β ELISA Kit | Catalog# EM0029 客户单位: 西湖大学

The cGAS-STING pathway activates transcription factor TFEB to stimulate lysosome biogenesis and pathogen clearance

期刊: Immunity IF: 25.5 DOI: 2025 Feb 11;58(2):309-325.e6.
产品名称: Human IL-6 ELISA Kit | Catalog# EH0001; Human IFN-β ELISA Kit | Catalog# EH0017; Human CXCL10 ELISA Kit | Catalog# EH0026
客户单位: 湖南第一师范大学

Inhalable Chemotactic Liposome for Targeted Modulation of Pulmonary Pre-Metastatic Niche

期刊: Adv. Funct. Mater. IF: 18.5 DOI: 2025 May 29;35(22):2419583.
产品名称: Human FGF basic ELISA Kit | Catalog# EH0024 客户单位: 四川大学

Biomimetic ferritin nanocages for synergistic co-delivery of metformin and rapamycin restore neurodevelopmental homeostasis in autism spectrum disorders

期刊: J Nanobiotechnology. IF: 12.6 DOI: 2025 Oct 14;23(1):670.
产品名称: Rat TNF-α ELISA Kit | Catalog# ER0007; Rat IL-6 ELISA Kit | Catalog# ER0004 客户单位: 上海大学

热销Quality ELISA 试剂盒清单精选

肿瘤免疫治疗			
Product name	Catalog#	Size	Sensitivity
Human alpha Fetoprotein ELISA Kit	EH0203	96T	6.88 pg/ml (50 µl); 32.49 pg/ml (10 µl)
Human B7-H2/CD275 ELISA Kit	EH0133	96T	0.01 pg/ml (50 µl); 4.95 pg/ml (10 µl)
Human B7-H6/NCR3LG1 ELISA Kit	EH0116	96T	3.98 pg/ml (50 µl); 4.61 pg/ml (10 µl)
Human BCMA/TNFRSF17 ELISA Kit	EH0143	96T	0.45 pg/ml (50 µl); 1.52 pg/ml (10 µl)
Human CDCP1 ELISA Kit	EH0083	96T	1.41 pg/ml (50 µl)
Human CTLA-4 ELISA Kit	EH0179	96T	1.33 pg/ml (50 µl); 1.97 pg/ml (10 µl)
Human EGF ELISA Kit	EH0023	96T	0.12 pg/ml (50 µl); 0.59 pg/ml (10 µl)
Human Endoglin/CD105 ELISA Kit	EH0204	96T	0.62 pg/ml (50 µl); 3.96 pg/ml (10 µl)
Human EphA2 ELISA Kit	EH0193	96T	0.30 pg/ml (50 µl); 1.23 pg/ml (10 µl)
Human FOLR1 ELISA Kit	EH0137	96T	2.40 pg/ml (50 µl); 3.97 pg/ml (10 µl)
Human G-CSF ELISA Kit	EH0149	96T	0.61 pg/ml (50 µl); 10.62 pg/ml (10 µl)
Human GM-CSF ELISA Kit	EH0051	96T	3.57 pg/ml (50 µl); 5.45 pg/ml (10 µl)
Human IDO1 ELISA Kit	EH0118	96T	1.61 pg/ml (50 µl); 3.04 pg/ml (10 µl)
Human M-CSF ELISA Kit	EH0126	96T	1.41 pg/ml (50 µl); 15.89 pg/ml (10 µl)
Human MERTK ELISA Kit (One step)	EHY0189	96T	0.03 ng/ml (50 µl); 0.41 ng/ml (10 µl)
Human Mesothelin ELISA Kit (One step)	EHY0152	96T	0.16 ng/ml (50 µl); 0.23 ng/ml (10 µl)
Human MICA ELISA Kit	EH0044	96T	8.3 pg/ml (50 µl)
Human MICB ELISA Kit	EH0127	96T	1.09 pg/ml (50 µl); 6.68 pg/ml (10 µl)
Human Osteoprotegerin/TNFRSF11B ELISA Kit	EH0070	96T	0.37 pg/ml (50 µl); 2.53 pg/ml (10 µl)
Human OX40/TNFRSF4 ELISA Kit	EH0120	96T	1.15 pg/ml (50 µl); 8.63 pg/ml (10 µl)
Human PD-1 ELISA Kit	EH0062	96T	0.05 pg/ml (50 µl); 1.64 pg/ml (10 µl)
Human PDGF-AA ELISA Kit	EH0058	96T	0.10 ng/ml (50 µl); 0.17 ng/ml (10 µl)
Human PDGF-BB ELISA Kit	EH0068	96T	2.79 pg/ml (50 µl); 7.56 pg/ml (10 µl)
Human TIM4 ELISA Kit (One step)	EHY0194	96T	0.001 ng/ml (50 µl); 0.03 ng/ml (10 µl)
Human VEGF ELISA Kit	EH0075	96T	1.04 pg/ml (50 µl); 3.21 pg/ml (10 µl)
Human VEGFR2/KDR ELISA Kit	EH0200	96T	1.54 pg/ml (50 µl); 9.53 pg/ml (10 µl)
Mouse GM-CSF/CSF2 ELISA Kit	EM0011	96T	0.19 pg/ml (10 µl)
Mouse PD-1 ELISA Kit	EM0037	96T	7.65 pg/ml (10 µl)
Mouse PD-L1 ELISA Kit	EM0036	96T	0.16 pg/ml (10 µl)
Mouse VEGF ELISA Kit	EM0013	96T	0.57 pg/ml (10 µl)

神经科学			
Product name	Catalog#	Size	Sensitivity
Human Alpha-Synuclein ELISA Kit	EH0191	96T	0.30 ng/ml (50 µl); 0.60 ng/ml (10 µl)
Human beta Amyloid 1-40 ELISA Kit	EH0111	96T	2.92 pg/ml (50 µl); 8.40 pg/ml (10 µl)
Human beta Amyloid 1-42 ELISA Kit	EH0112	96T	3.15 pg/ml (50 µl); 6.07 pg/ml (10 µl)
Human beta-NGF ELISA Kit	EH0065	96T	1.91 pg/ml (50 µl); 3.27 pg/ml (10 µl)
Human CD115/M-CSF R/CSF-1-R ELISA Kit (One step)	EHY0168	96T	0.01 ng/ml (50 µl); 0.08 ng/ml (10 µl)
Human CHIT1 ELISA Kit	EH0186	96T	0.22 pg/ml (50 µl); 0.49 pg/ml (10 µl)
Human Clusterin ELISA Kit (One step)	EHY0073	96T	0.14 ng/ml (10 µl)
Human CNTF ELISA Kit	EH0147	96T	0.25 pg/ml (50 µl); 3.77 pg/ml (10 µl)
Human E-Cadherin ELISA Kit	EH0115	96T	25.58 pg/ml (50 µl)
Human EphA2 ELISA Kit	EH0193	96T	0.30 pg/ml (50 µl); 1.23 pg/ml (10 µl)
Human FGF basic/FGF2/bFGF ELISA Kit	EH0024	96T	1.15 pg/ml (50 µl); 9.93 pg/ml (10 µl)
Human Lipocalin-2/NGAL ELISA Kit (One step)	EHY0166	96T	0.01 ng/ml (50 µl); 0.03 ng/ml (10 µl)
Human NCAM 1/CD56 ELISA Kit	EH0102	96T	3.69 pg/ml (50 µl); 37.50 pg/ml (10 µl)
Human Neurogranin ELISA Kit	EH0210	96T	0.92 pg/ml (50 µl); 11.03 pg/ml (10 µl)
Human Neuropilin-1/NRP1 ELISA Kit (One step)	EHY0175	96T	0.02 ng/ml (50 µl); 0.03 ng/ml (10 µl)
Human Oxytocin-neurophysin 1/OXT ELISA Kit	EH0153	96T	1.02 pg/ml (50 µl); 4.22 pg/ml (10 µl)
Human PDGF-BB ELISA Kit	EH0068	96T	2.79 pg/ml (50 µl); 7.56 pg/ml (10 µl)
Human Total HSP60 ELISA Kit	EH0060	96T	3.22 pg/ml (50 µl); 4.45 pg/ml (10 µl)
Human TREM2 ELISA Kit (One step)	EHY0122	96T	0.004 ng/ml (50 µl); 0.05 ng/ml (10 µl)
Human VCAM-1/CD106 ELISA Kit	EH0067	96T	7.07 pg/ml (50 µl)
Human VCAM-1/CD106 ELISA Kit (One step)	EHY0067	96T	0.04 ng/ml (50 µl); 0.08 ng/ml (10 µl)
Human YKL-40/CHI3L1 ELISA Kit (One step)	EHY0159	96T	0.01 ng/ml (50 µl); 0.05 ng/ml (10 µl)
Human/Mouse BDNF ELISA Kit	EH0069	96T	1.05 pg/ml (50 µl); 2.53 pg/ml (10 µl)
Mammal BACE-1 ELISA Kit	EH0195	96T	5.24 pg/ml (50 µl); 6.31 pg/ml (10 µl)
Mammal GFAP ELISA Kit	EH0178	96T	2.52 pg/ml (50 µl); 17.06 pg/ml (10 µl)
Mammal NF-L ELISA Kit	EH0154	96T	0.54 pg/ml (50 µl); 1.31 pg/ml (10 µl)
Mouse CRP/C-Reactive Protein ELISA Kit	EM0035	96T	0.01 ng/ml (50 µl)
Mouse TREM2 ELISA Kit	EM0039	96T	0.60 pg/ml (10 µl)
代谢性疾病			
Product name	Catalog#	Size	Sensitivity
Human Adiponectin ELISA Kit	EH0099	96T	0.15 pg/ml (50 µl); 3.56 pg/ml (10 µl)
Human Asprosin ELISA Kit	EH0032	96T	0.78 ng/ml (50 µl)

代谢性疾病			
Product name	Catalog#	Size	Sensitivity
Human CCL2/MCP-1 ELISA Kit	EH0010	96T	0.002pg/ml (50 µl); 0.02 pg/ml (10 µl)
Human FGF-19 ELISA Kit	EH0140	96T	0.48 pg/ml (50 µl); 4.11 pg/ml (10 µl)
Human Growth Hormone ELISA Kit	EH0142	96T	0.34 pg/ml (50 µl); 1.81 pg/ml (10 µl)
Human HGF ELISA Kit	EH0019	96T	4.07 pg/ml (50 µl) ; 6.01 pg/ml (10 µl)
Human IGFBP3/IBP3 ELISA Kit (One step)	EHY0157	96T	0.004 ng/ml (50 µl); 13.60 ng/ml (10 µl)
Human IL-1β ELISA Kit	EH0006	96T	0.09 pg/ml (50 µl) ; 0.18 pg/ml (10 µl)
Human IL-33 ELISA Kit	EH0047	96T	0.76 pg/ml (50 µl); 1.25 pg/ml (10 µl)
Human IL-4R ELISA Kit	EH0145	96T	0.24 pg/ml (50 µl); 2.90 pg/ml (10 µl)
Human IL-6 ELISA Kit	EH0001	96T	0.26 pg/ml (50 µl) ; 0.43 pg/ml (10 µl)
Human Leptin ELISA Kit	EH0056	96T	9.45 pg/ml (50 µl); 18.83 pg/ml (10 µl)
Human Leptin Receptor ELISA Kit	EH0187	96T	0.90 pg/ml (50 µl); 5.07 pg/ml (10 µl)
Human Lipocalin-2/NGAL ELISA Kit (One step)	EHY0166	96T	0.01 ng/ml (50µl); 0.03 ng/ml (10µl)
Human Lipoprotein A/LPA ELISA Kit	EH0171	96T	1.11 pg/ml (50 µl); 7.76 pg/ml (10 µl)
Human Prolactin ELISA Kit	EH0089	96T	3.47 pg/ml (50 µl); 11.20 pg/ml (10 µl)
Human Resistin ELISA Kit	EH0206	96T	0.33 pg/ml (50 µl); 0.48 pg/ml (10 µl)
Human SCF/KITLG ELISA Kit	EH0174	96T	0.22 pg/ml (50 µl); 2.44 pg/ml (10 µl)
Human TNF-α ELISA Kit	EH0002	96T	0.21 pg/ml (50 µl) ; 0.26 pg/ml (10 µl)
Human/Rat/Mouse Irisin ELISA Kit	EH0081	96T	2.60 pg/ml (50 µl); 9.39 pg/ml (10 µl)
Mammal Insulin ELISA Kit	EH0038	96T	60.07 pg/ml (50 µl) ; 314.56 pg/ml (10 µl)
Mammal VEGF-C ELISA Kit	EH0113	96T	1.82 pg/ml (50 µl)
Mouse IL-1β ELISA Kit	EM0029	96T	0.91 pg/ml (10 µl)
Mouse IL-1β Super Sensitivity ELISA Kit (高敏)	EMSS0029	96T	0.12 pg/ml (50 µl); 0.36 pg/ml (10 µl)
Mouse IL-33 ELISA Kit	EM0033	96T	0.33 pg/ml (10 µl)
Mouse IL-6 ELISA Kit	EM0004	96T	0.36 pg/ml (10 µl)
Mouse Leptin ELISA Kit	EM0027	96T	2.24 pg/ml (10 µl)
Mouse TNF-α ELISA Kit	EM0010	96T	0.097 pg/ml (10 µl)
Mouse VEGF ELISA Kit	EM0013	96T	0.57 pg/ml (10 µl)
Mouse/Rat IGF1 ELISA Kit	EMY0024	96T	0.10 pg/ml (50 µl); 0.37 pg/ml (10 µl)
Mouse/Rat IGF1/IGF-1/IGF-1 ELISA Kit	EM0024	96T	1.76 pg/ml (50 µl); 3.43 pg/ml (10 µl)
Rat IL-1β ELISA Kit	ER0001	96T	17.11 pg/ml (50 µl)
Rat IL-6 ELISA Kit	ER0004	96T	7.05pg/ml (50 µl)
Rat TNF-α ELISA Kit	ER0007	96T	0.08 pg/ml (50 µl); 0.16 pg/ml (10 µl)

心血管研究			
Product name	Catalog#	Size	Sensitivity
Human Angiogenin ELISA Kit (One step)	EHY0190	96T	0.01 ng/ml (50 µl); 0.07 ng/ml (10 µl)
Human CD31/PECAM-1 ELISA Kit	EH0192	96T	0.06 pg/ml (50 µl); 1.94 pg/ml (10 µl)
Human CD62E ELISA Kit	EH0198	96T	1.02 pg/ml (50 µl); 5.13 pg/ml (10 µl)
Human CKMB ELISA Kit	EH0072	96T	0.005 pg/ml (50 µl); 0.07 ng/ml (10 µl)
Human Complement C3 ELISA Kit (One step)	EHY0202	96T	0.09 pg/ml (10 µl)
Human CRP ELISA Kit	EH0076	96T	4.26 pg/ml (50 µl); 21.23 pg/ml (10 µl)
Human Endoglin/CD105 ELISA Kit	EH0204	96T	0.62 pg/ml (50 µl); 3.96 pg/ml (10 µl)
Human Growth Hormone ELISA Kit	EH0142	96T	0.34 pg/ml (50 µl); 1.81 pg/ml (10 µl)
Human HGF ELISA Kit	EH0019	96T	4.07 pg/ml (50 µl) ; 6.01 pg/ml (10 µl)
Human ICAM-1/CD54 ELISA Kit	EH0064	96T	0.49 pg/ml(50 µl); 6.36 pg/ml (10 µl)
Human ICAM-1/CD54 ELISA Kit (One step)	EHY0064	96T	1.22 pg/ml (50 µl); 0.05 ng/ml (10 µl)
Human Lipoprotein A/LPA ELISA Kit	EH0171	96T	1.11 pg/ml (50 µl); 7.76 pg/ml (10 µl)
Human PAI-1/Serpin E1 ELISA Kit (One step)	EHY0160	96T	2.18 pg/ml (50 µl); 40.00 pg/ml (10 µl)
Human PDGF-BB ELISA Kit	EH0068	96T	2.79 pg/ml (50 µl); 7.56 pg/ml (10 µl)
Human Proprotein Convertase 9/PCSK9 ELISA Kit	EH0071	96T	0.02 ng/ml (50 µl); 0.17 ng/ml (10 µl)
Human Proprotein Convertase 9/PCSK9 ELISA Kit (One step)	EHY0071	96T	0.39 ng/ml (50 µl); 0.73 ng/ml (10 µl)
Human RAGE ELISA Kit (One step)	EHY0222	96T	21.39 pg/ml (50 µl); 65.37 pg/ml (10 µl)
Human RBP4 ELISA Kit (One step)	EHY0212	96T	0.002 ng/ml (50 µl); 0.03 ng/ml (10 µl)
Human SAA1 ELISA Kit	EH0180	96T	0.30pg/ml (50 µl); 1.56 pg/ml (10 µl)
Human tPA ELISA Kit (One step)	EHY0158	96T	0.03 ng/ml (50 µl); 0.10 ng/ml (10 µl)
Human VCAM-1/CD106 ELISA Kit	EH0067	96T	7.07 pg/ml (50 µl)
Human VCAM-1/CD106 ELISA Kit (One step)	EHY0067	96T	0.04 ng/ml (50 µl); 0.08 ng/ml (10 µl)
Human VEGF ELISA Kit	EH0075	96T	1.04 pg/ml (50 µl); 3.21 pg/ml (10 µl)
Human VEGFR2/KDR ELISA Kit	EH0200	96T	1.54 pg/ml (50 µl); 9.53 pg/ml (10 µl)
Human von Willebrand Factor ELISA Kit (One step)	EHY0138	96T	0.02 ng/ml (50 µl); 0.69 ng/ml (10 µl)
Human/Mouse/Rat/Porcine/Canine TGF-β1 ELISA Kit	EH0012	96T	33.32 pg/ml (50 µl)
Mammal VEGF-C ELISA Kit	EH0113	96T	1.82 pg/ml (50 µl)
Mouse Coagulation Factor XII ELISA Kit	EMY0001	96T	0.22 ng/ml (10 µl)
Mouse Coagulation Factor XII ELISA Kit	EM0041	96T	2.35 pg/ml (10 µl)
Mouse CRP/C-Reactive Protein ELISA Kit	EM0035	96T	0.01 ng/ml (50 µl)
Mouse CXCL9/MIG ELISA Kit	EM0015	96T	0.29 pg/ml (10 µl)
Mouse VEGF ELISA Kit	EM0013	96T	0.57 pg/ml (10 µl)


心血管研究			
Product name	Catalog#	Size	Sensitivity
Rat CRP/C-Reactive Protein ELISA Kit	ER0009	96T	27.24 pg/ml (50 µl); 63.91 pg/ml (10 µl)
Rat CRP/C-Reactive Protein ELISA Kit (One step)	ERY0009	96T	2.9 pg/ml (50 µl); 42.15 pg/ml (10 µl)
骨科研究			
Product name	Catalog#	Size	Sensitivity
Human CD254/RANKL/TNFSF11 ELISA Kit	EH0214	96T	0.20 pg/ml (50 µl); 0.76 pg/ml (10 µl)
Human cIAP2 ELISA Kit	EH0123	96T	0.81 pg/ml (50 µl); 1.59 pg/ml (10 µl)
Human CTGF/CCN2 ELISA Kit (One step)	EHY0156	96T	0.07 ng/ml (50 µl); 0.35 ng/ml (10 µl)
Human IL-1β ELISA Kit	EH0006	96T	0.09 pg/ml (50 µl) ; 0.18 pg/ml (10 µl)
Human IL-6 ELISA Kit	EH0001	96T	0.26 pg/ml (50 µl) ; 0.43 pg/ml (10 µl)
Human Leptin ELISA Kit	EH0056	96T	9.45 pg/ml (50 µl); 18.83 pg/ml (10 µl)
Human Leptin Receptor ELISA Kit	EH0187	96T	0.90 pg/ml (50 µl); 5.07 pg/ml (10 µl)
Human Osteopontin (OPN) ELISA Kit	EH0057	96T	0.04 ng/ml (50 µl); 0.12ng/ml (10 µl)
Human Osteoprotegerin/TNFRSF11B ELISA Kit	EH0070	96T	0.37 pg/ml (50 µl); 2.53 pg/ml (10 µl)
Human TNF RII/TNFRSF1B ELISA Kit	EH0148	96T	0.27 pg/ml (50 µl); 5.77 pg/ml (10 µl)
Human TNF-α ELISA Kit	EH0002	96T	0.21 pg/ml (50 µl) ; 0.26 pg/ml (10 µl)
Human TRAF1 ELISA Kit	EH0215	96T	2.85 pg/ml (50 µl); 7.62 pg/ml (10 µl)
Mammal Insulin ELISA Kit	EH0038	96T	60.07 pg/ml (50 µl) ; 314.56 pg/ml (10 µl)
Mouse IL-1β ELISA Kit	EM0029	96T	0.91 pg/ml (10 µl)
Mouse IL-1β Super Sensitivity ELISA Kit (高敏)	EMSS0029	96T	0.12 pg/ml (50 µl); 0.36 pg/ml (10 µl)
Mouse IL-6 ELISA Kit	EM0004	96T	0.36 pg/ml (10 µl)
Mouse Leptin ELISA Kit	EM0027	96T	2.24 pg/ml (10 µl)
Mouse Pro-collagen I alpha 1/COL1A1 ELISA Kit	EMY0045	96T	0.01 pg/ml (10 µl)
Mouse RANKL/TRANCE/TNFSF11 ELISA Kit	EM0044	96T	0.36 pg/ml (10 µl)
Mouse TNF-α ELISA Kit	EM0010	96T	0.097 pg/ml (10 µl)
Rat IL-1β ELISA Kit	ER0001	96T	17.11 pg/ml (50 µl)
Rat IL-6 ELISA Kit	ER0004	96T	7.05pg/ml (50 µl)
Rat TNF-α ELISA Kit	ER0007	96T	0.08 pg/ml (50 µl); 0.16 pg/ml (10 µl)
肾脏病研究			
Product name	Catalog#	Size	Sensitivity
Human CCL2/MCP-1 ELISA Kit	EH0010	96T	0.002pg/ml (50 µl); 0.02 pg/ml (10 µl)
Human CD63 ELISA Kit	EH0167	96T	0.24 pg/ml (50 µl); 6.71 pg/ml (10 µl)
Human Clusterin ELISA Kit (One step)	EHY0073	96T	0.14 ng/ml (10 µl)

肾脏病研究			
Product name	Catalog#	Size	Sensitivity
Human CXCL1/GRO alpha ELISA Kit	EH0054	96T	14.89 pg/ml (50 µl); 32.67 pg/ml (10 µl)
Human EGF ELISA Kit	EH0023	96T	0.12 pg/ml (50 µl) ; 0.59 pg/ml (10 µl)
Human Endoglin/CD105 ELISA Kit	EH0204	96T	0.62 pg/ml (50 µl); 3.96 pg/ml (10 µl)
Human IL-10 ELISA Kit	EH0009	96T	0.89 pg/ml (10 µl)
Human IL-18 ELISA Kit	EH0053	96T	0.56pg/ml (50 µl); 1.28 pg/ml (10 µl)
Human IL-18BP ELISA Kit	EH0131	96T	1.67 pg/ml (50 µl); 1.70 pg/ml (10 µl)
Human IL-6 ELISA Kit	EH0001	96T	0.26 pg/ml (50 µl) ; 0.43 pg/ml (10 µl)
Human Lipocalin-2/NGAL ELISA Kit (One step)	EHY0166	96T	0.01 ng/ml (50µl); 0.03 ng/ml (10µl)
Human Osteoprotegerin/TNFRSF11B ELISA Kit	EH0070	96T	0.37 pg/ml (50 µl); 2.53 pg/ml (10 µl)
Human PDGF-AA ELISA Kit	EH0058	96T	0.10 ng/ml (50 µl); 0.17 ng/ml (10 µl)
Human TGF alpha ELISA Kit	EH0059	96T	1.91 pg/ml (50 µl); 3.89 pg/ml (10 µl)
Human TNF-α ELISA Kit	EH0002	96T	0.21 pg/ml (50 µl) ; 0.26 pg/ml (10 µl)
Human uPAR ELISA Kit	EH0077	96T	8.79 pg/ml (50 µl)
Human VCAM-1/CD106 ELISA Kit (One step)	EHY0067	96T	0.04 ng/ml (50 µl); 0.08 ng/ml (10 µl)
Human VEGF ELISA Kit	EH0075	96T	1.04 pg/ml (50 µl); 3.21 pg/ml (10 µl)
Mammal VEGF-C ELISA Kit	EH0113	96T	1.82 pg/ml (50 µl)
Mouse CCL2/MCP-1 ELISA Kit	EM0025	96T	1.78 pg/ml (10 µl)
Mouse CD68/SR-D1 ELISA Kit	EM0042	96T	0.62 pg/ml (10 µl)
Mouse CXCL1/KC ELISA Kit	EM0017	96T	2.11 pg/ml (10 µl)
Mouse IL-10 ELISA Kit	EM0005	96T	5.75 pg/ml (10 µl)
Mouse IL-18 ELISA Kit	EM0034	96T	1.46 pg/ml (10 µl)
Mouse IL-6 ELISA Kit	EM0004	96T	0.36 pg/ml (10 µl)
Mouse TNF-α ELISA Kit	EM0010	96T	0.097 pg/ml (10 µl)
Mouse VEGF ELISA Kit	EM0013	96T	0.57 pg/ml (10 µl)
Rat IL-10 ELISA Kit	ER0005	96T	5.76 pg/ml (50 µl)
Rat IL-6 ELISA Kit	ER0004	96T	7.05pg/ml (50 µl)
Rat TNF-α ELISA Kit	ER0007	96T	0.08 pg/ml (50 µl); 0.16 pg/ml (10 µl)
眼科研究			
Product name	Catalog#	Size	Sensitivity
Human CCL2/MCP-1 ELISA Kit	EH0010	96T	0.002pg/ml (50 µl); 0.02 pg/ml (10 µl)
Human Complement C3 ELISA Kit (One step)	EHY0202	96T	0.09 pg/ml (10 µl)
Human GM-CSF ELISA Kit	EH0051	96T	3.57 pg/ml (50 µl); 5.45 pg/ml (10 µl)

眼科研究

Product name	Catalog#	Size	Sensitivity
Human IFN- γ ELISA Kit	EH0039	96T	0.36 pg/ml (50 μ l); 0.85 pg/ml (10 μ l)
Human IL-6 ELISA Kit	EH0001	96T	0.26 pg/ml (50 μ l); 0.43 pg/ml (10 μ l)
Human MERTK ELISA Kit (One step)	EHY0189	96T	0.03 ng/ml (50 μ l); 0.41 ng/ml (10 μ l)
Human NG2/MCSP/CSPG4 ELISA Kit (One step)	EHY0162	96T	0.004 ng/ml (50 μ l); 0.21 ng/ml (10 μ l)
Human VE Cadherin ELISA Kit	EH0091	96T	0.01 ng/ml (50 μ l); 0.12 ng/ml (10 μ l)
Human VEGF ELISA Kit	EH0075	96T	1.04 pg/ml (50 μ l); 3.21 pg/ml (10 μ l)
Mammal VEGF-C ELISA Kit	EH0113	96T	1.82 pg/ml (50 μ l)
Mouse CCL2/MCP-1 ELISA Kit	EM0025	96T	1.78 pg/ml (10 μ l)
Mouse GM-CSF/CSF2 ELISA Kit	EM0011	96T	0.19 pg/ml (10 μ l)
Mouse IFN- γ ELISA Kit	EM0009	96T	3.24 pg/ml (10 μ l)
Mouse IL-6 ELISA Kit	EM0004	96T	0.36 pg/ml (10 μ l)
Mouse VEGF ELISA Kit	EM0013	96T	0.57 pg/ml (10 μ l)
Rat CXCL1/GRO alpha ELISA Kit	ER0010	96T	0.26 pg/ml (50 μ l); 0.46 pg/ml (10 μ l)
Rat IL-6 ELISA Kit	ER0004	96T	7.05pg/ml (50 μ l)

扫码了解更多
“Quality ELISA试剂盒”
产品信息

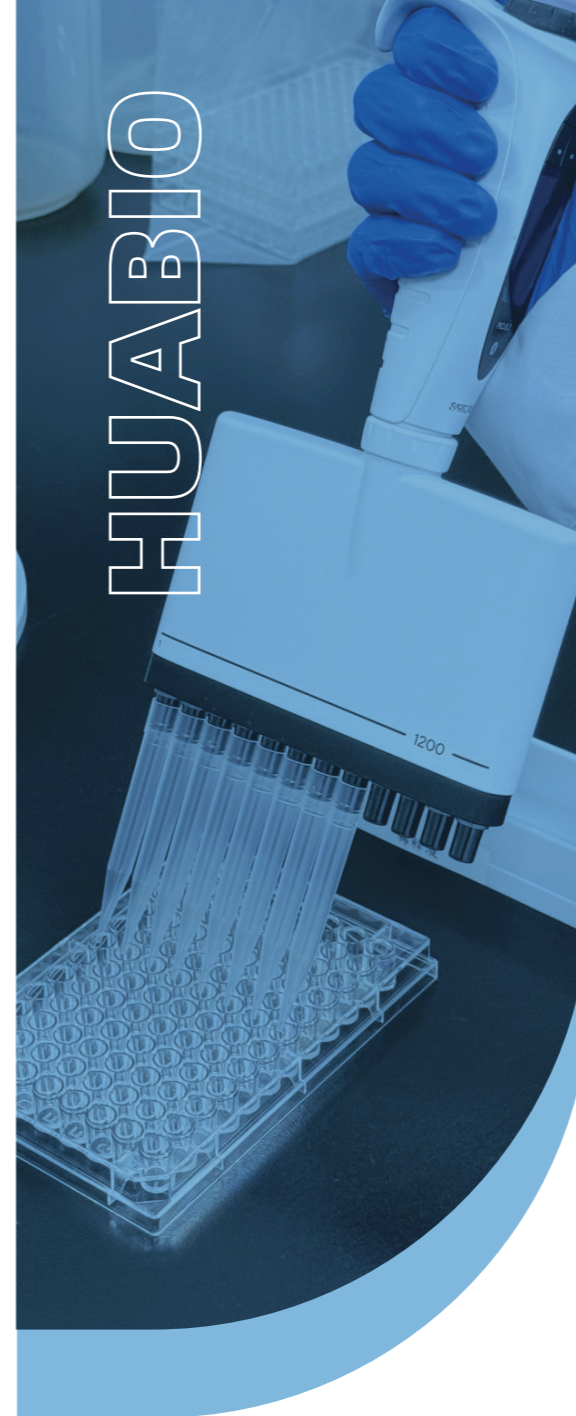

自主灵活： DTSet ELISA 试剂盒

华安生物 (HUABIO) 提供的DTSet ELISA 试剂盒包含夹心法免疫分析的必需组分，具有显著的自主灵活性和高性价比的特征。同时，能够兼顾应用场景的通用性及实验结果的稳定可靠。

试剂盒组分		Volume
捕获抗体	CA (capture antibody 100 \times)	1 vial
标准品	S (Standard)	1 vial
检测抗体	DA (Detect Antibody 100 \times)	1 vial
链霉亲和素-HRP	SH (Streptavidin-HRP 100 \times)	1 vial



产品优势

经深度验证的抗体对

每个试剂盒的核心抗体对均经过严格的配对筛选与性能验证，确保高特异性、低交叉反应性，并提供精准校准的标准品曲线

灵活的实验方案优化空间

支持根据样本特殊性自主优化孵育时间或稀释比例，以获取最优数据

适合大规模研究的高性价比

针对长期、多样本量的研究项目，该系列提供了更具成本效益的解决方案

支持定制化服务

提供针对新型或特定物种靶标的抗体对开发与试剂盒定制服务，助力前沿探索

热销DTSet ELISA 试剂盒清单精选

Product name	Catalog#	Size	Catalog#	Size	Sensitivity
Human IL-5 DTSet ELISA Kit	DSEH000805	96T*5	DSEH000815	96T*15	1.69 pg/ml(50 µl)
Human IL-2 DTSet ELISA Kit	DSEH001305	96T*5	DSEH001315	96T*15	0.25 pg/ml(50 µl)
Human IL-8 DTSet ELISA Kit	DSEH001405	96T*5	DSEH001415	96T*15	0.11 pg/ml(50 µl)
Human IFN-beta DTSet ELISA Kit	DSEH001705	96T*5	DSEH001715	96T*15	6.26 pg/ml(50 µl)
Human HGF DTSet ELISA Kit	DSEH001905	96T*5	DSEH001915	96T*15	0.45 pg/ml(50 µl)
Human MMP-9 DTSet ELISA Kit	DSEH002005	96T*5	DSEH002015	96T*15	1.27 pg/ml(50 µl)
Human EGF DTSet ELISA Kit	DSEH002305	96T*5	DSEH002315	96T*15	0.12 pg/ml(50 µl)
Human FGF-2 DTSet ELISA Kit	DSEH002405	96T*5	DSEH002415	96T*15	1.66 pg/ml(50 µl)
Human CXCL10/IP-10 DTSet ELISA Kit	DSEH002605	96T*5	DSEH002615	96T*15	0.008 pg/ml(50 µl)
Human IL-15 DTSet ELISA Kit	DSEH003105	96T*5	DSEH003115	96T*15	0.06 pg/ml(50 µl)
Human IFN-γ DTSet ELISA Kit	DSEH003905	96T*5	DSEH003915	96T*15	0.14 pg/ml(50 µl)
Human IL-7 DTSet ELISA Kit	DSEH004505	96T*5	DSEH004515	96T*15	0.86 pg/ml(50 µl)
Human CD85d DTSet ELISA Kit	DSEH004605	96T*5	DSEH004615	96T*15	77.03 pg/ml(50 µl)
Human IL-33 DTSet ELISA Kit	DSEH004705	96T*5	DSEH004715	96T*15	2.21 pg/ml(50 ul)
Human IL-13 DTSet ELISA Kit	DSEH004905	96T*5	DSEH004915	96T*15	10.26pg/ml(50 ul)
Human IL-3 DTSet ELISA Kit	DSEH005005	96T*5	DSEH005015	96T*15	0.02 pg/ml(50 µl)
Human GM-CSF DTSet ELISA Kit	DSEH005105	96T*5	DSEH005115	96T*15	0.06 pg/ml(50 µl)
Human IL-18 DTSet ELISA Kit	DSEH005305	96T*5	DSEH005315	96T*15	1.26 pg/ml(50 µl)
Human CXCL1/GRO alpha DTSet ELISA Kit	DSEH005405	96T*5	DSEH005415	96T*15	16.04 pg/ml(50 µl)
Human ICAM-1/CD54 DTSet ELISA Kit	DSEH006405	96T*5	DSEH006415	96T*15	2.48pg/ml(50 µl)
Human PDGF-BB DTSet ELISA Kit	DSEH006805	96T*5	DSEH006815	96T*15	0.10 pg/ml(50 µl)
Human CRP DTSet ELISA Kit	DSEH007605	96T*5	DSEH007615	96T*15	8.06 pg/ml(50 µl)
Human IFN-α DTSet ELISA Kit	DSEH008505	96T*5	DSEH008515	96T*15	0.06 pg/ml(50 µl)
Human IL-22 DTSet ELISA Kit	DSEH008605	96T*5	DSEH008615	96T*15	0.41pg/ml(50 µl)
Human CD200R1 DTSet ELISA Kit	DSEH008705	96T*5	DSEH008715	96T*15	0.24 pg/ml(50 µl)
Human Prolactin DTSet ELISA Kit	DSEH008905	96T*5	DSEH008915	96T*15	3.25 pg/ml(50 µl)
Human MBL2 DTSet ELISA Kit	DSEH009005	96T*5	DSEH009015	96T*15	7.53 pg/ml(50 µl)
Human VE Cadherin DTSet ELISA Kit	DSEH009105	96T*5	DSEH009115	96T*15	48.52pg/ml(50 µl)
Human IL-27 DTSet ELISA Kit	DSEH009405	96T*5	DSEH009415	96T*15	0.40 pg/ml(50 µl)
Human GITR/TNFRSF18 DTSet ELISA Kit	DSEH009505	96T*5	DSEH009515	96T*15	1.49 pg/ml(50 µl)
Human LIF DTSet ELISA Kit	DSEH009605	96T*5	DSEH009615	96T*15	0.15 pg/ml(50 µl)

Product name	Catalog#	Size	Catalog#	Size	Sensitivity
Human IL-31 DTSet ELISA Kit	DSEH009705	96T*5	DSEH009715	96T*15	0.07 pg/ml(50 µl)
Human IL-12B p40 DTSet ELISA Kit	DSEH009805	96T*5	DSEH009815	96T*15	0.03 pg/ml(50 µl)
Human Adiponectin DTSet ELISA Kit	DSEH009905	96T*5	DSEH009915	96T*15	0.23pg/ml(50 µl)
Human CD276 DTSet ELISA Kit	DSEH010105	96T*5	DSEH010115	96T*15	0.12pg/ml(50 µl)
Human MMP-8 DTSet ELISA Kit	DSEH010305	96T*5	DSEH010315	96T*15	3.33 pg/ml(50 µl)
Human IFN-γ R1 DTSet ELISA Kit	DSEH010405	96T*5	DSEH010415	96T*15	0.12pg/ml(50 µl)
Human CD30 DTSet ELISA Kit	DSEH010505	96T*5	DSEH010515	96T*15	0.74 pg/ml(50 µl)
Human CD33 DTSet ELISA Kit	DSEH010605	96T*5	DSEH010615	96T*15	0.17 pg/ml(50 µl)
Human HVEM/TNFRSF14 DTSet ELISA Kit	DSEH010705	96T*5	DSEH010715	96T*15	2.08pg/ml(50 µl)
Human Complement C5 DTSet ELISA Kit	DSEH010805	96T*5	DSEH010815	96T*15	0.80pg/ml(50 µl)
Human IL-23 p19 DTSet ELISA Kit	DSEH010905	96T*5	DSEH010915	96T*15	0.28pg/ml(50 µl)
Human beta Amyloid 1-40 DTSet ELISA Kit	DSEH011105	96T*5	DSEH011115	96T*15	0.55pg/ml(50 µl)
Human beta Amyloid 1-42 DTSet ELISA Kit	DSEH011205	96T*5	DSEH011215	96T*15	15.61 pg/ml(50 µl)
Mammal VEGF-C DTSet ELISA Kit	DSEH011305	96T*5	DSEH011315	96T*15	4.23 pg/ml(50 µl)
Human CD27/TNFRSF7 DTSet ELISA Kit	DSEH011405	96T*5	DSEH011415	96T*15	0.00pg/ml(50 µl)
Human E-Cadherin DTSet ELISA Kit	DSEH011505	96T*5	DSEH011515	96T*15	14.75pg/ml(50 µl)
Human B7-H6/NCR3LG1 ELISA Kit	DSEH011605	96T*5	DSEH011615	96T*15	13.41 pg/ml(50 µl)
Human CCL25 ELISA Kit	DSEH011705	96T*5	DSEH011715	96T*15	4.37 pg/ml(50 µl)
Human IDO1 DTSet ELISA Kit	DSEH011805	96T*5	DSEH011815	96T*15	3.78 pg/ml(50 µl)
Human CD48 DTSet ELISA Kit	DSEH011905	96T*5	DSEH011915	96T*15	0.67pg/ml(50 µl)
Human OX40/TNFRSF4 DTSet ELISA Kit	DSEH012005	96T*5	DSEH012015	96T*15	1.08 pg/ml(50 µl)
Human Human B7-1/CD80 DTSet ELISA Kit	DSEH012105	96T*5	DSEH012115	96T*15	0.22pg/ml(50 µl)
Human TREM2 DTSet ELISA Kit	DSEH012205	96T*5	DSEH012215	96T*15	0.46 pg/ml(50 µl)
Human cIAP2 DTSet ELISA Kit	DSEH012305	96T*5	DSEH012315	96T*15	0.11 pg/ml(50 µl)
Human CD137 DTSet ELISA Kit	DSEH012405	96T*5	DSEH012415	96T*15	1.55 pg/ml(10 µl)
Human MMP-7 DTSet ELISA Kit	DSEH012505	96T*5	DSEH012515	96T*15	10.53 pg/ml(50 µl)
Human M-CSF DTSet ELISA Kit	DSEH012605	96T*5	DSEH012615	96T*15	0.81 pg/ml(50 µl)
Human MICB DTSet ELISA Kit	DSEH012705	96T*5	DSEH012715	96T*15	6.98 pg/ml(50 µl)
Human TNF R1 DTSet ELISA Kit	DSEH012805	96T*5	DSEH012815	96T*15	0.00pg/ml(50 µl)
Human Granzyme A/GZMA DTSet ELISA Kit	DSEH012905	96T*5	DSEH012915	96T*15	0.08 pg/ml(50 µl)
Human CCL23/MPIF-1 ELISA Kit	DSEH013005	96T*5	DSEH013015	96T*15	0.06 pg/ml(50 µl)
Mouse IL-1 beta DTSet ELISA Kit	DSEM000105	96T*5	DSEM000115	96T*15	0.54pg/ml(10 µl)

Product name	Catalog#	Size	Catalog#	Size	Sensitivity
Mouse IL-4 DTSet ELISA Kit	DSEM000305	96T*5	DSEM000315	96T*15	0.81 pg/ml(10 µl)
Mouse IL-6 DTSet ELISA Kit	DSEM000405	96T*5	DSEM000415	96T*15	0.88pg/ml(10 µl)
Mouse IL-13 DTSet ELISA Kit	DSEM000705	96T*5	DSEM000715	96T*15	0.19 pg/ml(10 µl)
Mouse IL-17A DTSet ELISA Kit	DSEM000805	96T*5	DSEM000815	96T*15	0.35 pg/ml(10 µl)
Mouse TNF alpha DTSet ELISA Kit	DSEM001005	96T*5	DSEM001015	96T*15	1.10 pg/ml(10 µl)
Mouse IL-12/IL-23 p40 DTSet ELISA Kit	DSEM001205	96T*5	DSEM001215	96T*15	0.02pg/ml(10 µl)
Mouse VEGF DTSet ELISA Kit	DSEM001305	96T*5	DSEM001315	96T*15	3.76pg/ml(10 µl)
Mouse IL-9 DTSet ELISA Kit	DSEM002205	96T*5	DSEM002215	96T*15	0.02 pg/ml(10 µl)
Mouse IL-33 DTSet ELISA Kit	DSEM003305	96T*5	DSEM003315	96T*15	1.84pg/ml(10 µl)
Mouse PD-L1 DTSet ELISA Kit	DSEM003605	96T*5	DSEM003615	96T*15	0.03 pg/ml(10 µl)
Mouse PD1 DTSet ELISA Kit	DSEM003705	96T*5	DSEM003715	96T*15	7.33 pg/ml(10 µl)
Mouse LIF DTSet ELISA Kit	DSEM003805	96T*5	DSEM003815	96T*15	1.56 pg/ml(10 µl)
Mouse TREM2 DTSet ELISA Kit	DSEM003905	96T*5	DSEM003915	96T*15	0.10 pg/ml(10 µl)
Rat IL-1beta DTSet ELISA Kit	DSE000105	96T*5	DSE000115	96T*15	1.92 pg/ml(50 µl)
Rat IL-2 DTSet ELISA Kit	DSE000205	96T*5	DSE000215	96T*15	3.43 pg/ml(50 µl)
Rat IL-4 DTSet ELISA Kit	DSE000305	96T*5	DSE000315	96T*15	0.40 pg/ml(50 µl)
Rat IL-6 DTSet ELISA Kit	DSE000405	96T*5	DSE000415	96T*15	15.08 pg/ml(50 µl)
Rat IL-10 DTSet ELISA Kit	DSE000505	96T*5	DSE000515	96T*15	3.44 pg/ml(50 µl)
Rat IFN gamma DTSet ELISA Kit	DSE000605	96T*5	DSE000615	96T*15	6.05 pg/ml(50 µl)

扫码了解更多
“DTSet ELISA试剂盒”
产品信息



配对抗体

抗体对筛选“万里挑一”，品质保障，成就高质量的ELISA试剂盒

华安生物（HUABIO）严格遵循夹心法ELISA筛选和验证捕获抗体和检测抗体。为保证实际应用场景中检测结果的特异性和一致性，华安生物（HUABIO）使用生物样本进行筛选高品质和性能稳定的抗体对。同时，华安生物（HUABIO）提供的抗体具有高稳定性和批次间一致性，确保实验结果的稳定可靠和可重复。



神经退行性疾病

β-Amyloid 1-40 & 1-42, Tau & Phospho-Tau (T217), α-Synuclein, GFAP, TREM2, Neurogranin, BDNF, NF-L, S100B, AChE, c-Fos.....



自身免疫性疾病

IFN-α/γ, IL-6, IL-17, TNF-α, CXCL10, CCL2, BAFF, IL-23A, CTLA-4, ICOS, TSLP, CCL5, Complement C5, SAA.....



心血管疾病

ACE2, CRP, PTX3, ApoB, Lipoprotein A, VCAM-1, ICAM-1, MMP-9, TIMP1/2, PCSK9, CKMB, MPO, vWF.....



代谢性疾病

Insulin, IGF1/IGFBP, Leptin, Adiponectin, FGF-21, GDF15, Resistin, RBP4, ANGPTL3/4, PCSK9, DPP4, NLRP3, IL-1β.....



肿瘤研究

PD-L1/PD-1, CTLA-4, EGFR, VEGFA/VEGFR, AFP, CEACAM5, Mesothelin, CTGF, Osteopontin, TGF-β1.....



感染性疾病

SARS-CoV-2 Nucleocapsid, HIV1 p24, IL-6, CRP, IFN-γ, CXCL10, CD4/CD8, PD-1, Granzyme B, TNF-alpha, PTX-3, TREM-1, SAA1.....

300+ 靶点覆盖多种疾病领域

产品优势

生物样本验证

使用血清、血浆及细胞样本进行筛选和验证，保证实验结果的特异性和一致性

高灵敏度

对目标抗原的亲和力高，准确检测pg级别蛋白，适用于低表达天然蛋白的检测

高稳定性

在标准曲线范围内，严格控制CV值低于10%，保证实验结果的准确可靠

高批次间一致性

采用HiMab重组兔单抗研发生产，确保稳定持续供应

自行搭配灵活兼容


提供PBS-only形式，实现捕获和检测抗体自行搭配，随时偶联

生物样本筛选，保证实验结果的特异性和一致性

人/小鼠/大鼠 NF-L配对抗体

Capture antibody: HA723330; Detect antibody: HA723331


NF-L是一种分子量最小（68kDa）、含量最为丰富的神经丝亚单位，负责维持神经系统功能。NF-L通常作为多种神经系统疾病的生物标志物，在神经科学研究和临床诊断中具有重要价值。

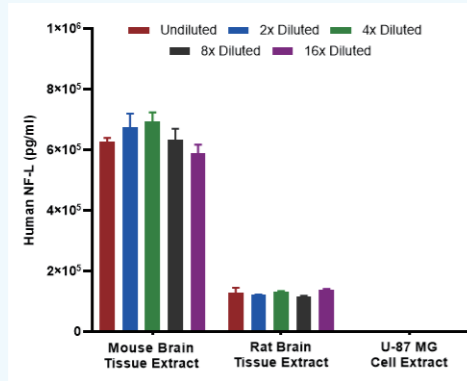

使用Human/Mouse/Rat NF-L配对抗体检测不同样本中NF-L的含量。结果表明：小鼠脑组织提取物上清中NF-L含量的平均值是643,742 pg/mL；大鼠脑组织提取物上清中NF-L含量的平均值是127,718 pg/m；U-87 MG细胞提取物上清中NF-L含量低于检测限。


检测灵敏度高，有效捕获pg/mL级别蛋白

Standard curve of human IL-3 matched pair antibodies

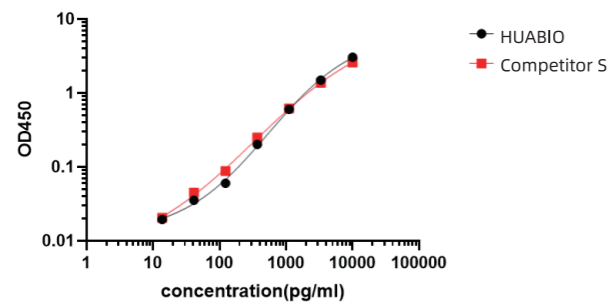



在相同实验操作条件下，使用Human IL-3配对抗体（Capture antibody: HA601153; Detect antibody: HA601155）和竞品S检测人源IL-3重组蛋白，并进行标准曲线统计分析。结果表明，较于竞品，华安生物（HUABIO）的IL-3配对抗体最低检测限为0.0140 pg/ml，检测灵敏度明显优于竞品。

高批次间一致性，保证实验数据的可重复性

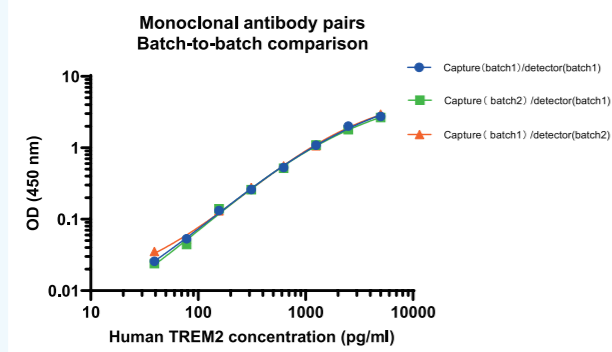


使用双抗夹心法对Human TREM2 配对抗体（Capture antibody: HA723149; Detect antibody: HA723151）的2次批次进行ELISA检测分析。结果表明，Human TREM2 配对抗体具有稳定的批次间一致性。

高稳定性，CV值低于10%，确保实验结果的稳定可靠

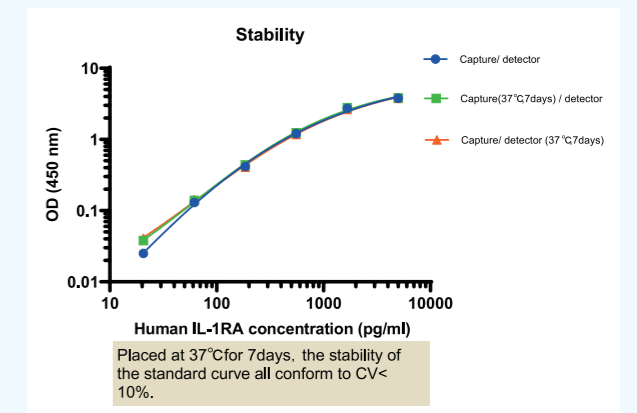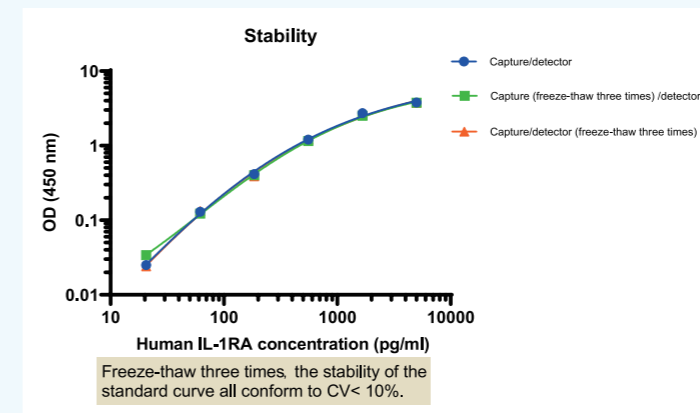

使用夹心ELISA方法对Human IL-1RA 配对抗体（Capture antibody: HA722788; Detect antibody: HA722789）分别对冻融稳定性（左图）和37°C热稳定性（右图）进行评价分析。结果表明，Human IL-1RA 配对抗体CV值严格控制在低于10%，具有高稳定性。

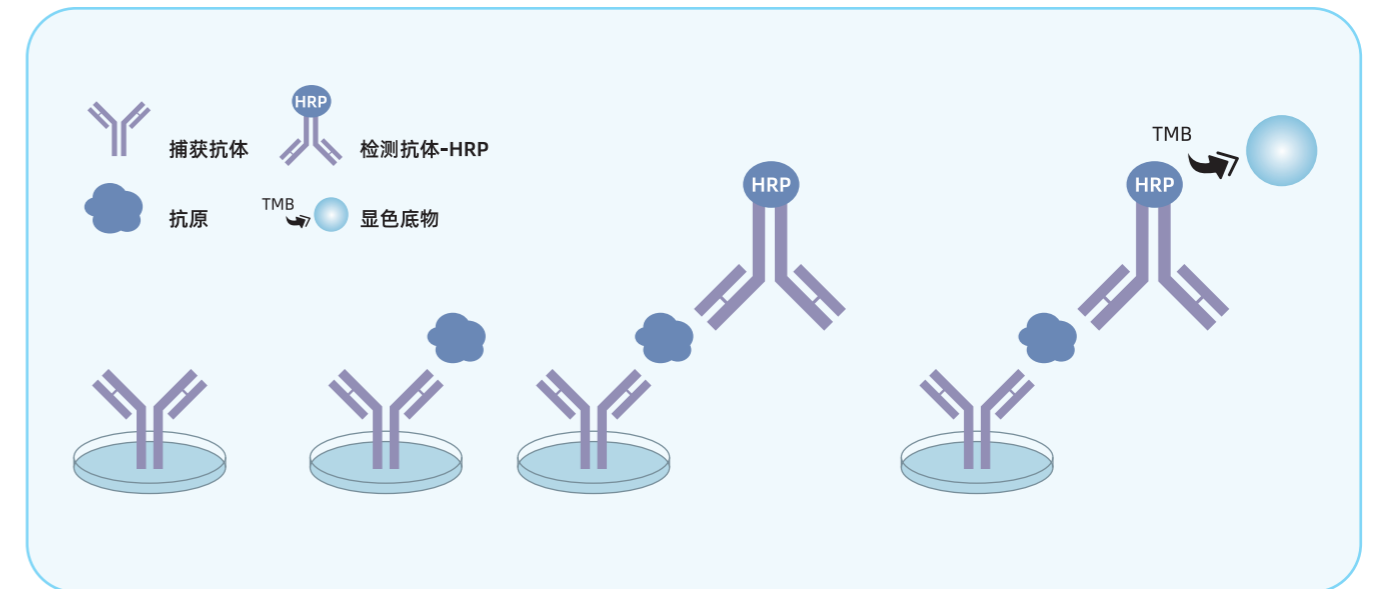

夹心法ELISA筛选配对抗体的实验流程

热销配对抗体清单精选

肿瘤免疫治疗				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
Human alpha Fetoprotein	HA723644	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human alpha Fetoprotein	HA723645	ELISA(Det)	H,Mk	Recombinant Rabbit mAb
Human B7-H6	HA723190	ELISA(Cap)	H	Recombinant Rabbit mAb
Human B7-H6	HA723191	ELISA(Det)	H	Recombinant Rabbit mAb
Human BCMA	HA725082	ELISA(Cap)	H	Recombinant Rabbit mAb
Human BCMA	HA725083	ELISA(Det)	H	Recombinant Rabbit mAb
Human CD147	HA725085	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD147	HA725086	ELISA(Det)	H	Recombinant Rabbit mAb
Human CD28	HA723457	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD28	HA723458	ELISA(Det)	H	Recombinant Rabbit mAb
Human CD30	HA722187	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD30	HA722188	ELISA(Det)	H	Recombinant Rabbit mAb
Human CD30-L	HA725339	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD30-L	HA725340	ELISA(Det)	H	Recombinant Rabbit mAb
Human CD33	HA722254	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD33	HA722255	ELISA(Det)	H	Recombinant Rabbit mAb
Human CD62L	HA723576	ELISA(Det),ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD62L	HA723578	ELISA(Det)	H	Recombinant Rabbit mAb
Human CD80	HA723123	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD80	HA723124	ELISA(Det)	H	Recombinant Rabbit mAb
Human CD85d	HA722381	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD85d	HA722382	ELISA(Det)	H	Recombinant Rabbit mAb
Human CEACAM5/CD66e	HA725349	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CEACAM5/CD66e	HA725350	ELISA(Det)	H	Recombinant Rabbit mAb
Human EGFR	HA724043	ELISA(Cap),ELISA(Det)	H	Recombinant Rabbit mAb
Human EGFR	HA724044	ELISA(Det)	H	Recombinant Rabbit mAb
Human ICOSLG/B7-H2/CD275	HA723320	ELISA(Cap)	H	Recombinant Rabbit mAb
Human ICOSLG/B7-H2/CD275	HA723321	ELISA(Det)	H	Recombinant Rabbit mAb
Human IDO1	HA723154	ELISA(Cap)	H	Recombinant Rabbit mAb
Human IDO1	HA723155	ELISA(Det)	H	Recombinant Rabbit mAb

肿瘤免疫治疗				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
Human Mesothelin	HA723287	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Mesothelin	HA723288	ELISA(Det)	H	Recombinant Rabbit mAb
Human Osteoprotegerin	HA725067	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Osteoprotegerin	HA725068	ELISA(Det)	H	Recombinant Rabbit mAb
Human PD1	HA723521	ELISA(Cap)	H	Recombinant Rabbit mAb
Human PD1	HA723522	ELISA(Det)	H	Recombinant Rabbit mAb
Human PDGF-AA	HA725061	ELISA(Cap)	H	Recombinant Rabbit mAb
Human PDGF-AA	HA725062	ELISA(Det)	H	Recombinant Rabbit mAb
Human PDGF-BB	HA722497	ELISA(Cap)	H	Recombinant Rabbit mAb
Human PDGF-BB	HA722498	ELISA(Det)	H	Recombinant Rabbit mAb
Human PD-L1	HA723740	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human PD-L1	HA723741	ELISA(Det)	H,Mk	Recombinant Rabbit mAb
Human PD-L2	HA723146	ELISA(Cap)	H	Recombinant Rabbit mAb
Human PD-L2	HA723147	ELISA(Det)	H	Recombinant Rabbit mAb
Human Poliovirus Receptor/PVR	HA723464	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Poliovirus Receptor/PVR	HA723465	ELISA(Det)	H	Recombinant Rabbit mAb
Human Prostate Specific Antigen	HA725257	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human Prostate Specific Antigen	HA725258	ELISA(Det)	H,Mk	Recombinant Rabbit mAb
Human Thymidine Kinase 1 / TK1	HA725239	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Thymidine Kinase 1 / TK1	HA725240	ELISA(Det)	H	Recombinant Rabbit mAb
Human TIM 4	HA723589	ELISA(Cap)	H	Recombinant Rabbit mAb
Human TIM 4	HA723590	ELISA(Det)	H	Recombinant Rabbit mAb
Human TIMP1	HA725367	ELISA(Cap)	H	Recombinant Rabbit mAb
Human TIMP1	HA725368	ELISA(Det)	H	Recombinant Rabbit mAb
Human VEGF Receptor 1	HA725365	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human VEGF Receptor 1	HA725366	ELISA(Det)	H,Mk	Recombinant Rabbit mAb
Human VEGFA	HA723493	ELISA(Cap)	H,D	Recombinant Rabbit mAb
Human VEGFA	HA723494	ELISA(Det)	H,D	Recombinant Rabbit mAb
Human/Mouse VEGFC	HA723241	ELISA(Cap)	H,Ms	Recombinant Rabbit mAb
Human/Mouse VEGFC	HA723242	ELISA(Det)	H,Ms	Recombinant Rabbit mAb
Mouse Osteopontin	HA723924	ELISA(Cap)	Ms	Recombinant Rabbit mAb
Mouse Osteopontin	HA723925	ELISA(Det)	Ms	Recombinant Rabbit mAb

肿瘤免疫治疗				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
Mouse PD1	HA722100	ELISA(Cap)	Ms	Recombinant Rabbit mAb
Mouse PD1	HA722101	ELISA(Det)	Ms	Recombinant Rabbit mAb
Mouse PD-L1	HA722189	ELISA(Cap)	Ms	Recombinant Rabbit mAb
Mouse PD-L1	HA722190	ELISA(Det)	Ms	Recombinant Rabbit mAb
Mouse VEGF	HA721850	ELISA(Cap)	Ms	Recombinant Rabbit mAb
Mouse VEGF	HA721851	ELISA(Det)	Ms	Recombinant Rabbit mAb
Mouse/Rat VEGFR3 / FLT4	HA724110	ELISA(Cap)	Ms,R	Recombinant Rabbit mAb
Mouse/Rat VEGFR3 / FLT4	HA724111	ELISA(Det)	Ms,R	Recombinant Rabbit mAb
TGF beta 1	HA723725	ELISA(Cap)	H,Ms,R,D	Recombinant Rabbit mAb
TGF beta 1	HA723726	ELISA(Det)	H,Ms,R,D	Recombinant Rabbit mAb
神经科学				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
beta Amyloid 1-40	HA722077	ELISA(Cap)	H,Ms	Recombinant Rabbit mAb
beta Amyloid 1-40	HA722078	ELISA(Det)	H,Ms	Recombinant Rabbit mAb
beta Amyloid 1-42	HA722079	ELISA(Cap)	H	Recombinant Rabbit mAb
beta Amyloid 1-42	HA722080	ELISA(Det)	H	Recombinant Rabbit mAb
Brain-derived tau/BD-tau	HA723847	ELISA(Cap)	H,Ms,R	Recombinant Rabbit mAb
Brain-derived tau/BD-tau	HA723848	ELISA(Det)	H	Recombinant Rabbit mAb
GFAP	HA601342	ELISA(Det)	H,Ms,R	Recombinant Mouse mAb
GFAP	HA722234	ELISA(Cap)	H,Ms,R	Recombinant Rabbit mAb
Human Acetylcholinesterase	HA725223	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Acetylcholinesterase	HA725224	ELISA(Det)	H	Recombinant Rabbit mAb
Human Alpha-Synuclein	HA723467	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Alpha-Synuclein	HA723468	ELISA(Det)	H	Recombinant Rabbit mAb
Human BACE1	HA725171	ELISA(Cap)	H	Recombinant Rabbit mAb
Human BACE1	HA725172	ELISA(Det)	H	Recombinant Rabbit mAb
Human Beta Nerve Growth Factor	HA722903	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Beta Nerve Growth Factor	HA722904	ELISA(Det)	H	Recombinant Rabbit mAb
Human CNTF	HA723381	ELISA(Cap)	H,R	Recombinant Rabbit mAb
Human CNTF	HA723382	ELISA(Det)	H,R	Recombinant Rabbit mAb
Human/Mouse c-Fos	HA723957	ELISA(Cap)	H,Ms	Recombinant Rabbit mAb
Human/Mouse c-Fos	HA723958	ELISA(Det)	H,Ms	Recombinant Rabbit mAb

神经科学				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
Human FAM171A2	HA724137	ELISA(Cap)	H	Recombinant Rabbit mAb
Human FAM171A2	HA724138	ELISA(Det)	H	Recombinant Rabbit mAb
Human LAMC2	HA723717	ELISA(Cap)	H	Recombinant Rabbit mAb
Human LAMC2	HA723718	ELISA(Det)	H	Recombinant Rabbit mAb
Human NCAM1	HA722946	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human NCAM1	HA722947	ELISA(Det),ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human Neurogranin	HA725191	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Neurogranin	HA725192	ELISA(Det)	H	Recombinant Rabbit mAb
Human Neuropilin-1	HA725112	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Neuropilin-1	HA725113	ELISA(Det)	H	Recombinant Rabbit mAb
Human NPTXR	HA725189	ELISA(Cap)	H	Recombinant Rabbit mAb
Human NPTXR	HA725190	ELISA(Det)	H	Recombinant Rabbit mAb
Human NPTX2	HA725187	ELISA(Cap)	H	Recombinant Rabbit mAb
Human NPTX2	HA725188	ELISA(Det)	H	Recombinant Rabbit mAb
Human Oxytocin-neurophysin 1	HA723452	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Oxytocin-neurophysin 1	HA723453	ELISA(Det)	H	Recombinant Rabbit mAb
Human PARK7	HA620001	ELISA(Cap)	H	Recombinant Mouse mAb
Human PARK7	HA620002	ELISA(Det)	H	Recombinant Mouse mAb
Human TREM2	HA723150	ELISA(Cap)	H	Recombinant Rabbit mAb
Human TREM2	HA723151	ELISA(Det)	H	Recombinant Rabbit mAb
Mouse/Rat Galectin 1	HA724222	ELISA(Cap)	Ms,R	Recombinant Rabbit mAb
Mouse/Rat Galectin 1	HA724223	ELISA(Det)	Ms,R	Recombinant Rabbit mAb
Human/Mouse/Rat NF-L	HA723330	ELISA(Cap)	H,Ms,R	Recombinant Rabbit mAb
Human/Mouse/Rat NF-L	HA723331	ELISA(Det)	H,Ms,R	Recombinant Rabbit mAb
Human REST	HA725193	ELISA(Cap)	H	Recombinant Rabbit mAb
Human REST	HA725194	ELISA(Det)	H	Recombinant Rabbit mAb
Human/Mouse/Rat Synapsin I	HA723993	ELISA(Cap)	H,Ms,R	Recombinant Rabbit mAb
Human/Mouse/Rat Synapsin I	HA723994	ELISA(Det)	H,Ms,R	Recombinant Rabbit mAb
Human Apolipoprotein E	HA724129	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Apolipoprotein E	HA724130	ELISA(Det)	H	Recombinant Rabbit mAb
Phospho-Tau (T217)	HA723842	ELISA(Cap)	H	Recombinant Rabbit mAb
Tau	HA723898	ELISA(Det)	H	Recombinant Rabbit mAb

神经科学				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
Phospho-Tau (T217)	HA723841	ELISA(Cap)	H,Ms	Recombinant Rabbit mAb
Phospho-Tau (T217)	HA723877	ELISA(Cap)	H	Recombinant Rabbit mAb
Phospho-Tau (T217)	HA723922	ELISA(Cap)	H,Ms	Recombinant Rabbit mAb
Tau	HA723209	ELISA(Det)	H,Ms,R	Recombinant Rabbit mAb
MTBR-TAU 243/cleaved tau (Val256)	HA724134	ELISA(Cap)	H	Recombinant Rabbit mAb
MTBR-TAU 243/cleaved tau (Val256)	HA724135	ELISA(Det)	H	Recombinant Rabbit mAb
S100B	HA723728	ELISA(Cap)	H,Ms,R	Recombinant Rabbit mAb
S100B	HA723729	ELISA(Det)	H,Ms,R	Recombinant Rabbit mAb
SNAP25	HA725227	ELISA(Cap)	H,Ms	Recombinant Rabbit mAb
SNAP25	HA725228	ELISA(Det)	H,Ms	Recombinant Rabbit mAb
Tau	HA723850	ELISA(Cap)	H,Ms,R	Recombinant Rabbit mAb
Tau	HA723851	ELISA(Det)	H	Recombinant Rabbit mAb
Human YKL-40	HA725088	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human YKL-40	HA725089	ELISA(Det)	H,Mk	Recombinant Rabbit mAb
心血管研究				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
H-FABP	HA725177	ELISA(Cap)	H,Ms,Mk	Recombinant Rabbit mAb
H-FABP	HA725178	ELISA(Det)	H,Ms,Mk	Recombinant Rabbit mAb
Human Angiogenin	HA725076	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Angiogenin	HA725077	ELISA(Det)	H	Recombinant Rabbit mAb
Human ANGPTL3	HA725079	ELISA(Cap)	H	Recombinant Rabbit mAb
Human ANGPTL3	HA725080	ELISA(Det)	H	Recombinant Rabbit mAb
Human Antithrombin III / ATIII	HA723612	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Antithrombin III / ATIII	HA723613	ELISA(Det)	H	Recombinant Rabbit mAb
Human Apolipoprotein A1	HA724026	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Apolipoprotein A1	HA724027	ELISA(Det)	H	Recombinant Rabbit mAb
Human CD31	HA722906	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD31	HA722907	ELISA(Det)	H	Recombinant Rabbit mAb
Human CKMB	HA723348	ELISA(Cap)	H,R,Ms,Mk	Recombinant Rabbit mAb
Human CKMB	HA723349	ELISA(Det)	H,R,Ms,Mk	Recombinant Rabbit mAb
Human Fetuin A	HA724046	ELISA(Cap),ELISA(Det)	H	Recombinant Rabbit mAb
Human Fetuin A	HA724047	ELISA(Det)	H	Recombinant Rabbit mAb

心血管研究				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
Human MMP13	HA723293	ELISA(Cap)	H	Recombinant Rabbit mAb
Human MMP13	HA723294	ELISA(Det)	H	Recombinant Rabbit mAb
Human MMP-3	HA722877	ELISA(Cap)	H	Recombinant Rabbit mAb
Human MMP-3	HA722878	ELISA(Det)	H	Recombinant Rabbit mAb
Human MMP-7	HA722926	ELISA(Cap)	H	Recombinant Rabbit mAb
Human MMP-7	HA722927	ELISA(Det)	H	Recombinant Rabbit mAb
Human MMP-8	HA721978	ELISA(Cap)	H	Recombinant Rabbit mAb
Human MMP-8	HA721979	ELISA(Det)	H	Recombinant Rabbit mAb
Human MMP-9	HA722106	ELISA(Cap)	H	Recombinant Rabbit mAb
Human MMP-9	HA722107	ELISA(Det)	H	Recombinant Rabbit mAb
Human MMP-9	HA722162	ELISA(Cap)	H	Recombinant Rabbit mAb
Human MMP-9	HA722163	ELISA(Det)	H	Recombinant Rabbit mAb
Human Myeloperoxidase	HA723673	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Myeloperoxidase	HA723674	ELISA(Det)	H	Recombinant Rabbit mAb
Human PAI1	HA723435	ELISA(Cap),ELISA(Det)	H	Recombinant Rabbit mAb
Human PAI1	HA723436	ELISA(Det)	H	Recombinant Rabbit mAb
Human PCSK9	HA725124	ELISA(Cap)	H	Recombinant Rabbit mAb
Human PCSK9	HA725125	ELISA(Det)	H	Recombinant Rabbit mAb
Human PDGF-BB	HA722497	ELISA(Cap)	H	Recombinant Rabbit mAb
Human PDGF-BB	HA722498	ELISA(Det)	H	Recombinant Rabbit mAb
Human Plasminogen	HA724132	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Plasminogen	HA724133	ELISA(Det)	H	Recombinant Rabbit mAb
Human RBP4	HA723688	ELISA(Cap)	H	Recombinant Rabbit mAb
Human RBP4	HA723689	ELISA(Det)	H	Recombinant Rabbit mAb
Human Renin	HA725169	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Renin	HA725170	ELISA(Det)	H	Recombinant Rabbit mAb
Human TIMP1	HA725367	ELISA(Cap)	H	Recombinant Rabbit mAb
Human TIMP1	HA725368	ELISA(Det)	H	Recombinant Rabbit mAb
Human VCAM-1	HA725064	ELISA(Cap)	H	Recombinant Rabbit mAb
Human VCAM-1	HA725065	ELISA(Det)	H	Recombinant Rabbit mAb
Human VEGF Receptor 1	HA725365	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human VEGF Receptor 1	HA725366	ELISA(Det)	H,Mk	Recombinant Rabbit mAb

心血管研究				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
Human VEGF Receptor 2	HA725121	ELISA(Cap)	H	Recombinant Rabbit mAb
Human VEGF Receptor 2	HA725122	ELISA(Det)	H	Recombinant Rabbit mAb
Human VEGFA	HA723493	ELISA(Cap)	H,D	Recombinant Rabbit mAb
Human VEGFA	HA723494	ELISA(Det)	H,D	Recombinant Rabbit mAb
Human Von Willebrand Factor	HA723058	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Von Willebrand Factor	HA723059	ELISA(Det)	H	Recombinant Rabbit mAb
Human/Monkey TIMP4	HA724220	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human/Monkey TIMP4	HA724221	ELISA(Det)	H,Mk	Recombinant Rabbit mAb
Human/Mouse VEGFC	HA723241	ELISA(Cap)	H,Ms	Recombinant Rabbit mAb
Human/Mouse VEGFC	HA723242	ELISA(Det)	H,Ms	Recombinant Rabbit mAb
Human/Mouse/Rat FGF1	HA723987	ELISA(Cap)	H,Ms,R,Mk	Recombinant Rabbit mAb
Human/Mouse/Rat FGF1	HA723988	ELISA(Det)	H,Ms,R,Mk	Recombinant Rabbit mAb
Mouse Angiotensinogen	HA725131	ELISA(Cap)	Ms	Recombinant Rabbit mAb
Mouse Angiotensinogen	HA725132	ELISA(Det)	Ms	Recombinant Rabbit mAb
Mouse Apolipoprotein B/Apo B	HA723944	ELISA(Cap)	Ms	Recombinant Rabbit mAb
Mouse Apolipoprotein B/Apo B	HA723945	ELISA(Det)	Ms	Recombinant Rabbit mAb
Mouse Apolipoprotein E	HA725221	ELISA(Cap)	Ms	Recombinant Rabbit mAb
Mouse Apolipoprotein E	HA725222	ELISA(Det)	Ms	Recombinant Rabbit mAb
代谢性疾病				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
AGRP	HA725307	ELISA(Cap)	H,Ms	Recombinant Rabbit mAb
AGRP	HA725308	ELISA(Det)	H,Ms	Recombinant Rabbit mAb
Human Adiponectin	HA722260	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human Adiponectin	HA722261	ELISA(Det)	H,Mk	Recombinant Rabbit mAb
Human ANGPTL3	HA725079	ELISA(Cap)	H	Recombinant Rabbit mAb
Human ANGPTL3	HA725080	ELISA(Det)	H	Recombinant Rabbit mAb
Human ANGPTL4	HA725389	ELISA(Cap)	H	Recombinant Rabbit mAb
Human ANGPTL4	HA725390	ELISA(Det)	H	Recombinant Rabbit mAb
Human DPP4 / CD26	HA724058	ELISA(Cap)	H	Recombinant Rabbit mAb
Human DPP4 / CD26	HA724059	ELISA(Det)	H	Recombinant Rabbit mAb
Human Fetuin A	HA724046	ELISA(Cap),ELISA(Det)	H	Recombinant Rabbit mAb
Human Fetuin A	HA724047	ELISA(Det)	H	Recombinant Rabbit mAb

代谢性疾病				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
Human FGF-19	HA725100	ELISA(Cap)	H	Recombinant Rabbit mAb
Human FGF-19	HA725101	ELISA(Det)	H	Recombinant Rabbit mAb
Human FGL1	HA724032	ELISA(Cap)	H	Recombinant Rabbit mAb
Human FGL1	HA724033	ELISA(Det)	H	Recombinant Rabbit mAb
Human GDF15	HA723805	ELISA(Cap)	H	Recombinant Rabbit mAb
Human GDF15	HA723806	ELISA(Det)	H	Recombinant Rabbit mAb
Human Growth Hormone	HA725106	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Growth Hormone	HA725107	ELISA(Det)	H	Recombinant Rabbit mAb
Human HGF	HA721778	ELISA(Cap)	H	Recombinant Rabbit mAb
Human HGF	HA721779	ELISA(Det)	H	Recombinant Rabbit mAb
Human Leptin	HA723421	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Leptin	HA723422	ELISA(Det)	H	Recombinant Rabbit mAb
Human Leptin Receptor	HA725058	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Leptin Receptor	HA725059	ELISA(Det)	H	Recombinant Rabbit mAb
Human Lipocalin-2	HA725109	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human Lipocalin-2	HA725110	ELISA(Det)	H,Mk	Recombinant Rabbit mAb
Human Lipoprotein A	HA723570	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Lipoprotein A	HA723571	ELISA(Det)	H	Recombinant Rabbit mAb
Human Prealbumin	HA725259	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Prealbumin	HA725260	ELISA(Det)	H	Recombinant Rabbit mAb
Human Prolactin	HA721941	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Prolactin	HA721942	ELISA(Det)	H	Recombinant Rabbit mAb
Human RAGE	HA725289	ELISA(Cap)	H	Recombinant Rabbit mAb
Human RAGE	HA725290	ELISA(Det)	H	Recombinant Rabbit mAb
Human RBP4	HA723688	ELISA(Cap)	H	Recombinant Rabbit mAb
Human RBP4	HA723689	ELISA(Det)	H	Recombinant Rabbit mAb
Human Resistin	HA723658	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human Resistin	HA723659	ELISA(Det)	H,Mk	Recombinant Rabbit mAb
Insulin	HA722096	ELISA(Cap)	H,Ms,R	Recombinant Rabbit mAb
Insulin	HA722097	ELISA(Det)	H,Ms,R	Recombinant Rabbit mAb
Mouse Adiponectin	HA723679	ELISA(Cap)	Ms	Recombinant Rabbit mAb
Mouse Adiponectin	HA723680	ELISA(Det)	Ms	Recombinant Rabbit mAb

代谢性疾病				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
Mouse FGF-21	HA723731	ELISA(Cap)	Ms	Recombinant Rabbit mAb
Mouse FGF-21	HA723732	ELISA(Det)	Ms	Recombinant Rabbit mAb
Mouse Lipocalin-2	HA725399	ELISA(Cap)	Ms	Recombinant Rabbit mAb
Mouse Lipocalin-2	HA725400	ELISA(Det)	Ms	Recombinant Rabbit mAb
Mouse Prolactin	HA723955	ELISA(Cap)	Ms	Recombinant Rabbit mAb
Mouse Prolactin	HA723956	ELISA(Det)	Ms	Recombinant Rabbit mAb
Mouse Resistin	HA723808	ELISA(Cap)	Ms	Recombinant Rabbit mAb
Mouse Resistin	HA723809	ELISA(Det)	Ms	Recombinant Rabbit mAb
Mouse/Rat GDF15	HA723866	ELISA(Cap)	Ms,R	Recombinant Rabbit mAb
Mouse/Rat GDF15	HA723867	ELISA(Det)	Ms,R	Recombinant Rabbit mAb
Mouse/Rat Leptin	HA723046	ELISA(Cap)	Ms,R	Recombinant Rabbit mAb
Mouse/Rat Leptin	HA723047	ELISA(Det)	Ms,R	Recombinant Rabbit mAb
感染/自身免疫				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
HIV1 p24	HA721527	ELISA(Cap)	Hiv1	Recombinant Rabbit mAb
HIV1 p24	HA721528	ELISA(Det)	Hiv1	Recombinant Rabbit mAb
Human CCL11/Eotaxin	HA725197	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human CCL11/Eotaxin	HA725198	ELISA(Det)	H,Mk	Recombinant Rabbit mAb
Human CCL16	HA725281	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CCL16	HA725282	ELISA(Det)	H	Recombinant Rabbit mAb
Human CCL18	HA724087	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CCL18	HA724088	ELISA(Det)	H	Recombinant Rabbit mAb
Human CD200	HA722379	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD200	HA722380	ELISA(Det)	H	Recombinant Rabbit mAb
Human CD200R1	HA722256	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD200R1	HA722257	ELISA(Det)	H	Recombinant Rabbit mAb
Human CD27/TNFRSF7	HA722973	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD27/TNFRSF7	HA722974	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD276	HA722491	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD276	HA722492	ELISA(Det)	H	Recombinant Rabbit mAb
Human CD30	HA722187	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CD30	HA722188	ELISA(Det)	H	Recombinant Rabbit mAb

感染/自身免疫				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
Human Complement C3	HA725140	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Complement C3	HA725141	ELISA(Det)	H	Recombinant Rabbit mAb
Human Complement C5	HA722682	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Complement C5	HA722683	ELISA(Det)	H	Recombinant Rabbit mAb
Human Complement factor B	HA725161	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Complement factor B	HA725162	ELISA(Det)	H	Recombinant Rabbit mAb
Human Complement factor D	HA724090	ELISA(Cap),ELISA(Det)	H	Recombinant Rabbit mAb
Human Complement factor D	HA724091	ELISA(Det)	H	Recombinant Rabbit mAb
Human CRP	HA722108	ELISA(Cap)	H,Mk,D	Recombinant Rabbit mAb
Human CRP	HA722109	ELISA(Det)	H,Mk,D	Recombinant Rabbit mAb
Human CXCL10 / IP-10	HA725007	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human CXCL10 / IP-10	HA725008	ELISA(Det)	H,Mk	Recombinant Rabbit mAb
Human CXCL13	HA722791	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CXCL13	HA722792	ELISA(Det)	H	Recombinant Rabbit mAb
Human CXCL13	HA722794	ELISA(Det)	H	Recombinant Rabbit mAb
Human CXCL16	HA725363	ELISA(Cap)	H	Recombinant Rabbit mAb
Human CXCL16	HA725364	ELISA(Det)	H	Recombinant Rabbit mAb
Human DR3/TNFRSF25	HA723869	ELISA(Cap)	H	Recombinant Rabbit mAb
Human DR3/TNFRSF25	HA723870	ELISA(Det)	H	Recombinant Rabbit mAb
Human Fas / CD95	HA723835	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Fas / CD95	HA723836	ELISA(Det)	H	Recombinant Rabbit mAb
Human Fas Ligand	HA724023	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Fas Ligand	HA724024	ELISA(Det)	H	Recombinant Rabbit mAb
Human FGL1	HA724032	ELISA(Cap)	H	Recombinant Rabbit mAb
Human FGL1	HA724033	ELISA(Det)	H	Recombinant Rabbit mAb
Human Flt3 Ligand / FLT3LG	HA725325	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Flt3 Ligand / FLT3LG	HA725326	ELISA(Det)	H	Recombinant Rabbit mAb
Human Galectin 3	HA725385	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Galectin 3	HA725386	ELISA(Det)	H	Recombinant Rabbit mAb
Human Galectin 9	HA725333	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Galectin 9	HA725334	ELISA(Det)	H	Recombinant Rabbit mAb
Human GITR	HA722383	ELISA(Cap)	H	Recombinant Rabbit mAb

感染/自身免疫				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
Human GITR	HA722384	ELISA(Det)	H	Recombinant Rabbit mAb
Human GM-CSF	HA721288	ELISA(Cap)	H,Mk,D	Recombinant Rabbit mAb
Human GM-CSF	HA721289	ELISA(Det)	H,Mk,D	Recombinant Rabbit mAb
Human Granzyme K	HA724061	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Granzyme K	HA724062	ELISA(Det)	H	Recombinant Rabbit mAb
Human IL-1 beta	HA721278	ELISA(Cap)	H,D	Recombinant Rabbit mAb
Human IL-1 beta	HA721279	ELISA(Det)	H,D	Recombinant Rabbit mAb
Human IL-10	HA721236	ELISA(Cap)	H,Mk,D	Recombinant Rabbit mAb
Human IL-10	HA721237	ELISA(Det)	H,Mk,D	Recombinant Rabbit mAb
Human IL-13	HA725013	ELISA(Cap)	H	Recombinant Rabbit mAb
Human IL-13	HA725014	ELISA(Det)	H	Recombinant Rabbit mAb
Human IL-17A	HA721232	ELISA(Cap)	H,Mk,D	Recombinant Rabbit mAb
Human IL-17A	HA721233	ELISA(Det)	H,Mk,D	Recombinant Rabbit mAb
Human IL-4	HA721230	ELISA(Cap)	H,Mk	Recombinant Rabbit mAb
Human IL-4	HA721231	ELISA(Det)	H,Mk	Recombinant Rabbit mAb
Human Interferon gamma	HA721684	ELISA(Cap)	H,D	Recombinant Rabbit mAb
Human Interferon gamma	HA721685	ELISA(Det)	H,D,Mk	Recombinant Rabbit mAb
Human NCR1/NKp46	HA725321	ELISA(Cap)	H	Recombinant Rabbit mAb
Human NCR1/NKp46	HA725322	ELISA(Det)	H	Recombinant Rabbit mAb
Human NKG2D	HA725299	ELISA(Cap)	H	Recombinant Rabbit mAb
Human NKG2D	HA725300	ELISA(Det)	H	Recombinant Rabbit mAb
Human NLRP3	HA723275	ELISA(Cap)	H	Recombinant Rabbit mAb
Human NLRP3	HA723276	ELISA(Det)	H	Recombinant Rabbit mAb
Human Pentraxin 3 / PTX3	HA723387	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Pentraxin 3 / PTX3	HA723388	ELISA(Det)	H	Recombinant Rabbit mAb
Human S100A12	HA725215	ELISA(Cap)	H	Recombinant Rabbit mAb
Human S100A12	HA725216	ELISA(Det)	H	Recombinant Rabbit mAb
Human SAA1	HA601333	ELISA(Cap)	H	Mouse mAb
Human SAA1	HA601334	ELISA(Det)	H	Mouse mAb
Human Siglec 7	HA723392	ELISA(Cap)	H	Recombinant Rabbit mAb
Human Siglec 7	HA723393	ELISA(Cap),ELISA(Det)	H	Recombinant Rabbit mAb
Human Siglec 9	HA723397	ELISA(Cap)	H	Recombinant Rabbit mAb

感染/自身免疫				
Product name	Catalog#	Tested applications	Reactivity	Antibody type
Human Siglec 9	HA723398	ELISA(Det)	H	Recombinant Rabbit mAb
Human TRAF1	HA723618	ELISA(Cap)	H	Recombinant Rabbit mAb
Human TRAF1	HA723619	ELISA(Det)	H	Recombinant Rabbit mAb
Human TREM1	HA725233	ELISA(Cap)	H	Recombinant Rabbit mAb
Human TREM1	HA725234	ELISA(Det)	H	Recombinant Rabbit mAb
Human TSLP	HA723487	ELISA(Cap)	H	Recombinant Rabbit mAb
Human TSLP	HA723488	ELISA(Det)	H	Recombinant Rabbit mAb
Human/Monkey Albumin	HA724127	ELISA(Cap)	H,Mk,D	Recombinant Rabbit mAb
Human/Monkey Albumin	HA724128	ELISA(Det)	H,Mk,D	Recombinant Rabbit mAb
Human/Mouse MASP2	HA725391	ELISA(Cap)	H,Ms	Recombinant Rabbit mAb
Human/Mouse MASP2	HA725392	ELISA(Det)	H,Ms	Recombinant Rabbit mAb
SARS-CoV-2 Nucleocapsid protein	HA600007	WB,ELISA(Cap)	Sars-cov-2	Mouse mAb
SARS-CoV-2 Nucleocapsid protein	HA600008	WB,ELISA(Cap)	Sars-cov-2	Mouse mAb

*H: Human; Ms: Mouse; R: Rat; Mk: Monkey; D: Dog

扫码了解更多
“配对抗体研究”
产品信息

